



ANALIZA BIAŁEK



100012878

## Laboratory of Analysis

**Client:** EMPEROR  
**Order number:** 100012878  
**Name:** EMPEROR TURINABOL  
**Number of samples:** 1 container

**Batch number:** -----  
**Recive date:** 2026-02-26  
**Start date:** 2026-03-13  
**End date:** 2026-03-16

## Sample identification

Substance	Method of taking a sample	Form
4-Chlorodehydromethyltestosterone	Sample provided by customer	White tablets in a plastic container


## Results

Substance	Sample name	Content	Purity
4-Chlorodehydromethyltestosterone	EMPEROR TURINABOL 10MG/TAB	10,14mg/tab	-----

## Research comment

The test sample was applied in the amount of 3 $\mu$ l to the BioBasic 8 Thermo Scientific column. During the course of the analysis with the method a linear gradient in the solutions of H<sub>2</sub>O and Acetonitrile with the addition of 0.1% trifluoroacetic acid, a peak of the sought substance appeared. Its concentration, when compared to the reference chromatogram, is 10,14mg/tab. Height the measurement error is + - 0,05 mg. The actual state of the sample sent before the test is presented on the photo below

## Photo of the test sample

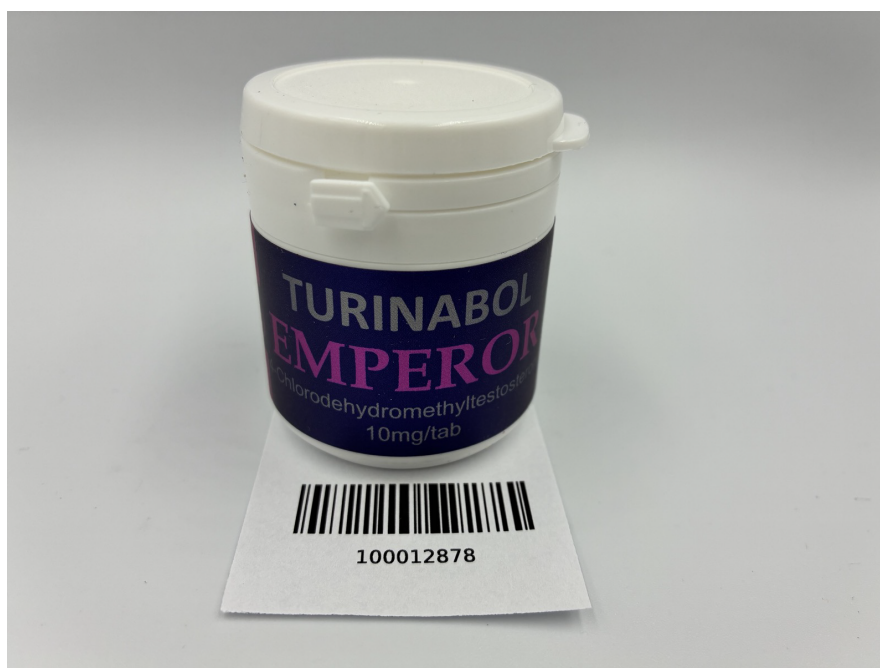

## Password

5SL9386H


**Analiza Białek sp. z o.o.**  
**KRS:** 0000966168  
**NIP:** 8992922970  
**REGON:** 521751875

**Headquarters and Laboratory:**  
ul. Duńska 9, 54-427 Wrocław  
Budynek DELTA

**Lab technician:**  
Analysis and research specialist  
mgr Arkadiusz Zajęc



## Chromatogram

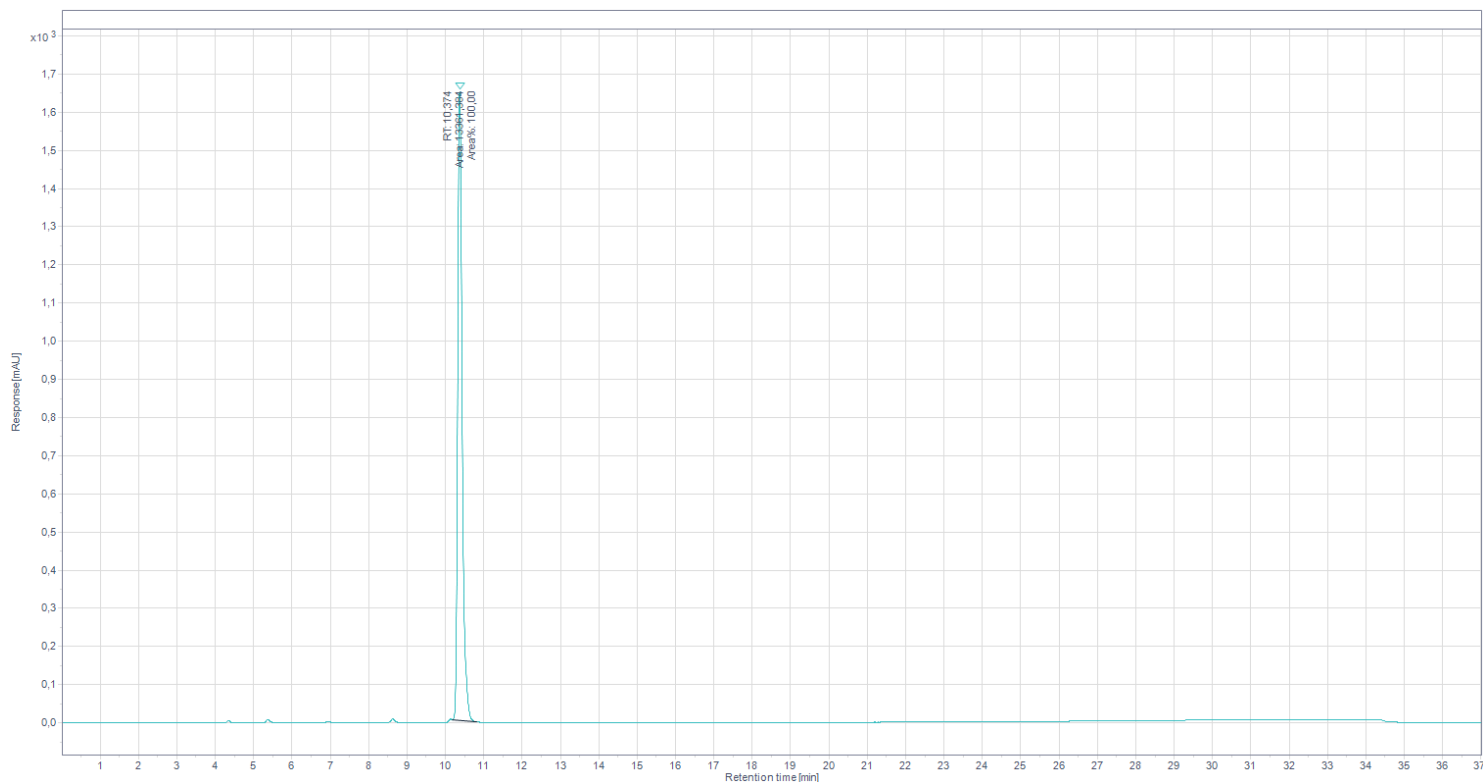

## Method

Time	H2O + 0,1% TFA	ACN + 0,1% TFA
0min	60%	40%
2min	60%	40%
27min	5%	95%
30min	5%	95%

## Password

5SL9386H